
©2008 Human Touch®, LLC. Patents pending and intellectual property rights enforced.
No medical claims are warranted or implied by the use of this product.
Human Touch is a trademark of Human Touch®, LLC.
HT-1260-NA-A0

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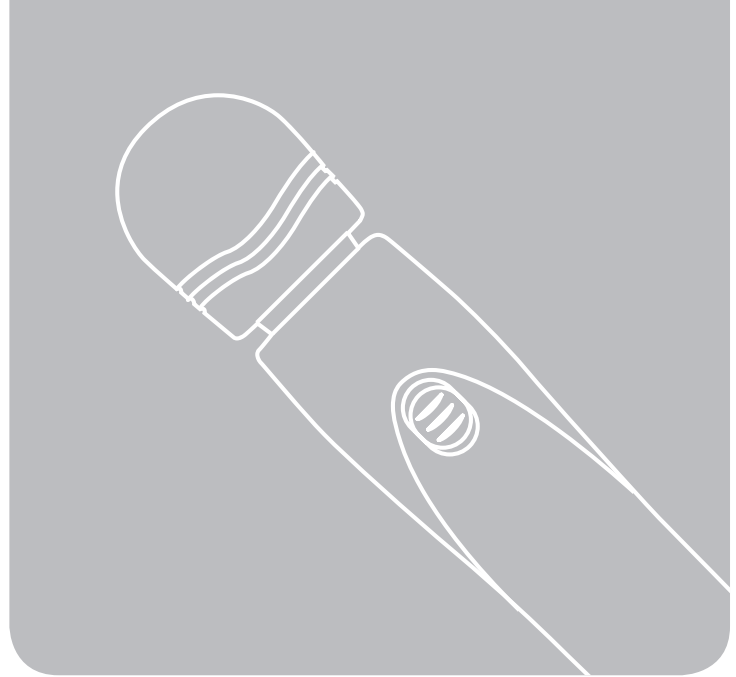ht¹²⁶⁰

acuvibe®
portable personal massager

Your Acuvibe® makes it easy to feel better whenever and wherever you'd like – its compact, cordless design makes it perfectly portable. Relax . . . invigorate . . . loosen muscles before exercise, or soothe stiffness afterward. However you use it, your Acuvibe is a perfect way to rejuvenate. Before using your Acuvibe, be sure to read this manual carefully.

keep it safe

IMPORTANT SAFETY INSTRUCTIONS

When using an electrical appliance, basic precautions should always be followed, including the following:

Read all instructions before using this appliance.

 **DANGER: To reduce the risk of electrical shock:**

- Always unplug this product from the electrical outlet immediately after using and before cleaning.
- Do not reach for this product if it has fallen into water. Unplug immediately.
- Do not use while bathing or in the shower.
- Do not place or store appliance where it can fall or be pulled into a tub or sink. Do not place into or drop into water or other liquid.

 **WARNING: To reduce the risk of burns, fire, electric shock or injury to persons:**

- Do not use on sensitive skin areas or in the presence of poor circulation.
- An appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
- Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.

- Close supervision is necessary when this product is used by, on, or near children, invalids, or disabled persons.
- Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return this product to a service center for examination and repair.
- Do not carry this appliance by supply cord or use cord as a handle.
- Always keep the power cord away from heated surfaces.
- Never operate this product with the air openings blocked. Keep the air openings free of lint, hair, and the like. Never operate on a soft surface such as a bed or couch where the air openings may be blocked.
- Never drop or insert any object into any opening.
- While using the massager, keep your hair away from the air inlets.
- Do not leave unattended on any surface while it is operating.
- Do not use outdoors.
- Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- To disconnect, turn all controls to the OFF position, then remove plug from outlet.
- Do not massage abdomen when pain is present or if you are pregnant.
- Do not massage any area of the body that is swollen, inflamed or covered with skin eruptions.
- Do not use on unexplained pain.
- Do not use this product when standing on a damp floor, or when any part of the body is in contact with plumbing or any similar ground.

- Do not use this product if noise above normal sound is heard.
- Never use while sleeping.
- Always avoid contact with the eye or eyelid.
- Never recharge unit using any charger other than the charger that came with the unit.
- Never expose batteries to excessive heat.

Safety Precautions

- Comfortable clothing is recommended while using this product.
- If you feel any abnormal discomfort during the operation of this product, turn off the power immediately and do not use the product. Consult your physician.
- When using this product after it has been in storage, check to see that the product is operating properly prior to use.

SAVE THESE INSTRUCTIONS

Do not attempt to repair this product yourself.

This product is designed for personal in-home use. Commercial use will void warranty.

Contains NiCd (Nickel-Cadmium) battery, must recycle or dispose of properly. Never discard any appliance containing batteries in the trash, as they are environmentally toxic when disposed of in landfills. Please contact your local recycling center for proper disposal information.

Care and Maintenance

If your massager becomes soiled, unplug the unit and wipe it clean with a cloth dampened with a non-abrasive cleanser. Do not use thinner, benzene or any other solvents. The product is specially designed and constructed to be maintenance free. No periodic lubrication or servicing is required. Your product should only be serviced by an authorized service center. Care should be taken not to overload the unit.

STORAGE:

Your product has been carefully designed to give you many years of trouble-free service. Always unplug the unit when not in use. Additionally, if not in use for a long period of time, we recommend that the power cord be coiled and that the unit be placed in an environment free of dust and moisture. Do not store near heat or open flame. Do not leave this product exposed to direct sunlight for long periods of time. This may cause fading or damage.

get ready



IMPORTANT: Charge your Acuvibe for eight hours before using it the first time.

1. Make sure the power switch is set to OFF.
2. Connect the charger cord to the charging jack at the end of the Acuvibe, then plug the charger into a wall outlet. The charging LED lights, indicating that the Acuvibe is charging.

NOTE: The total run time is dependent upon the charge level, temperature (lower temperatures result in shorter run times) and the amount of pressure applied during the massage.

IMPORTANT: Completely drain the battery and re-charge it for 16 hours once each month.

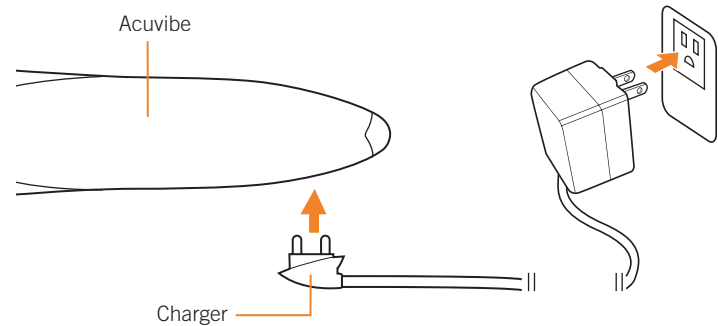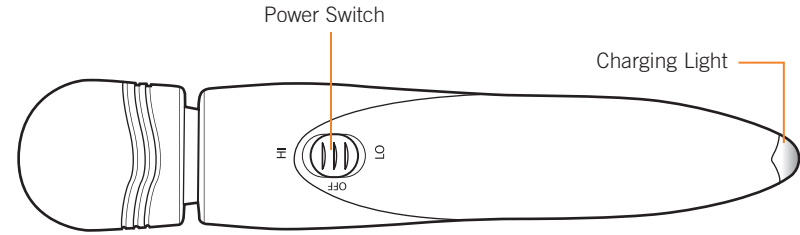take control

1. Disconnect the charger from your Acuvibe then turn the power switch to LO or HI.

NOTE: *The massager will not work while plugged into the charger.*

2. Use the soft, spring-mounted "flex-node" massage head to soothe muscle stiffness away. For extra massage benefit, gently stretch the muscle as you are massaging it.

CAUTION: It's best to use the massager on bare skin. Do not massage any single area for more than 15 minutes. Do not massage the upper neck area, near the base of the skull.


3. When you have completed your massage, turn the power switch to OFF.

 *If the Acuvibe overheats, it will shut down automatically and remain inoperable until it is safe to use again.*

Specifications*

Transformer voltage:	Input 120V 60Hz 25W Output 3 VDC 1500mA
Massager voltage:	2.4 VDC 15W
Transformer weight:	1.1 Lb.
Massager weight:	1.2 Lb.
Material:	ABS and EVA

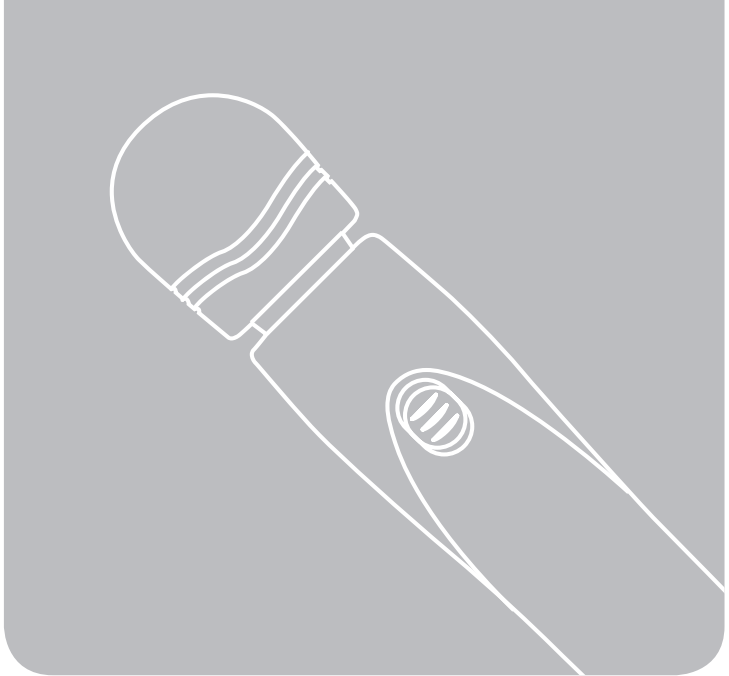**Specifications are subject to change without prior notice.*

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