




Violight Personal/Travel Toothbrush Sanitizer.
A toothbrush case with built in UV power. Free Toothbrush included! \$29.95

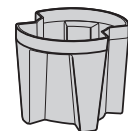

Replacement Drip Cup (Set of 2). Keep your countertop Violight extra fresh with a new drip cup! \$9.95

Violight Dental Spa. UV + Sonic cleaning for dentures, retainers, mouthguards & more! \$59.95



Violight UV Wand. Disinfect any surface with a wave of the wand. Great for countertops, keyboards, phones, and toys! \$79.95

Get FREE Ground Shipping at www.violight.com! Enter Coupon Code "SHIP" at checkout!

SANTIZING YOUR TOOTHBRUSH CAN BE SO EASY!

Congratulations on your purchase of the VIOlight Toothbrush Sanitizer and Storage System. As a health-conscious consumer, you're aware of the benefits of keeping a clean bathroom environment, changing your toothbrush every 3 months, and the importance of sanitizing it every day.

With your new VIOlight Toothbrush Sanitizer and Storage System, sanitizing your toothbrush is so easy. There are virtually no extra steps to take beyond your normal brushing routine. You brush, rinse, insert your toothbrush into the VIOlight, press the button and you're sanitizing. See the blue glow, and you'll know it's working! Within just 10 minutes, the germicidal ultraviolet light will eliminate up to 99% of the germs it comes into contact with.

The removable drip cup allows for easy cleaning and its four separate slots prevent your toothbrushes from touching and cross-contaminating each other. The VIOpeg inserts let you sanitize electric toothbrush heads, children's toothbrushes, or any toothbrush with short handles.

Best of all, your toothbrushes are stored in an enclosed, sanitized environment, clean and ready for the next time you brush.

By following the instructions below, you and your family will enjoy safe, sanitized toothbrushes for many years to come.

Joel Pinsky
Chairman and CEO

Jonathan Pinsky
President



toothbrush sanitizer & storage system

Model VIO100



SAFEGUARDS AND PRECAUTIONS IMPORTANT TO READ CAREFULLY

- As with any electrical appliance, care should be taken in its use and proper instruction and supervision should be given to children.
- Do not place the unit too close to sinks, toilets or bathtubs where it can be in danger of falling in.
- If unit falls into water, do not attempt to remove until you unplug from wall socket.
- Do not immerse the base or fill it with liquid of any kind.
- Always unplug the unit when cleaning or not in use.
- Do not use unit if the cord is frayed or the plug damaged.
- Keep cord away from hot surfaces.
- The unit is equipped with safety mechanisms that will not allow the bulb to operate if the cap is off.
- If ever the bulb stays on when the cap is off, never look directly at the lit bulb. Immediately unplug and return for service.
- Do not use your unit for any other purpose other than sanitizing toothbrushes.
- If unit is dropped or damaged or if the bulb remains on after the 10 minute sanitizing cycle, return it to the place where purchased or directly to VIOlight for inspection or repair.
- Do not insert toothbrushes or operate the unit without the drip cup being in place.
- Only use the power adapter that comes with the unit.





INSTRUCTIONS FOR USE

1. Remove all components from the box.
2. To open the VIOlight, hold the base in one hand and with your other hand, gently rotate the frosted cap CLOCKWISE to release the safety lock.
3. Remove the cap.
4. Carefully remove any protective material and the bulb from the inside of the unit. Using a tissue or soft cloth, screw the bulb securely into its socket at the bottom center of the unit by gently turning clockwise. ***As soon as you feel resistance, STOP turning. Do not over-tighten the UV bulb.***
5. Lift the removable drip cup straight up and out of the base, and rinse with warm water before first use. Re-insert the drip cup to its original position, firmly seated in the base.
6. If you will be using **electric toothbrush heads**, use the VIOpegs. Simply insert the VIOpegs into the drip cup slots that will be utilized for the electric toothbrush heads.
7. Replace the frosted cap onto the top of the base unit. Line up the cap's safety guides with the slots on top of the base and gently turn the cap COUNTERCLOCKWISE until you hear a secure click.
8. Plug the power supply cord into the base unit, and insert the two-prong plug into your 110V wall outlet.

The VIOlight is now ready for sanitizing.
Before inserting toothbrushes, please shake off excess water from the bristles.

9. Carefully insert your toothbrush (only one per slot) into the holes on into the VIOlight frosted cap. Be sure the bristles face the center of the unit. If you are using electric toothbrush heads, be sure to place them into the slot(s) that contain the VIOpegs. (See step 6 above)
10. On the base of the unit press the blue VIO button once, and the UV light will turn on and the frosted cap will glow blue. Your toothbrushes are now sanitizing. The light will automatically shut off after 10 minutes, leaving your toothbrushes fully sanitized. If it is ever necessary to turn off the unit before sanitization is complete, simply press the blue VIO button again.

The VIOlight has a built in safety lock feature that will not allow operation any time the cap is off or if the cap is not safely rotated and securely locked in place.

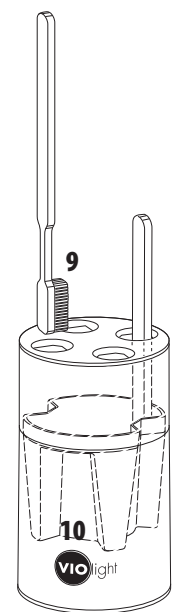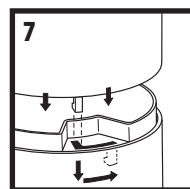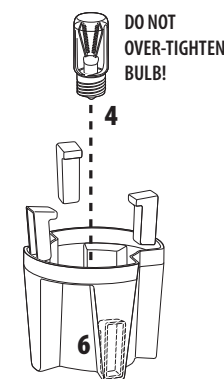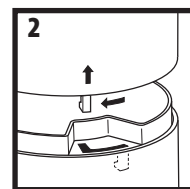

CLEANING & CARE OF YOUR VIOLIGHT

To reduce the risk of electrical shock or injury always unplug the unit before cleaning or replacing the bulb.

We recommend you rinse the removable drip cup at least once a week. More frequent rinsing may be necessary depending on usage. Please check your drip cup for residue!

Remove all toothbrushes from the unit. Follow step 3 and step 5 above to remove and clean the drip cup.

If necessary, use a moist cloth or sponge to clean the VIOlight base as needed. NEVER immerse the base in water, or place under running water for any reason. Do not fill the base with any liquid.

The frosted cap and the VIOpegs can also be rinsed in warm water whenever necessary.

REPLACING THE UV BULB

To reduce the risk of electrical shock or injury always unplug the unit before cleaning or replacing the bulb.

The UV Germicidal bulb has a life span of over 2,000 hours. As with any bulb, efficiency decreases with time. We recommended you change the bulb at least once a year.

Remove the drip cup and using a soft cloth or tissue unscrew the bulb by gently rotating counterclockwise.

Do not squeeze the bulb too hard or bend from side to side. **Carefully review step 4 when inserting your new UV bulb.**

Bulbs and other replacement parts can be ordered through www.violight.com or by calling 888-882-2367, Mon-Fri 9am-5pm.

TROUBLESHOOTING

If for any reason your unit fails to work after pressing the on/off button, check for the following:

- Did you press the VIO button hard enough?
 - Is the transformer cord securely attached at both ends?
 - Unplug and check if the bulb is loose and needs tightening or if the filament is broken and bulb needs to be replaced.
 - Is cap on securely and properly locked into position?
 - Does a wall outlet switch have to be turned on?
- Call Violight customer service at 888-882-2367 Mon-Fri 9am-5pm or visit www.violight.com

ONE YEAR LIMITED WARRANTY

Violight Inc. warrants to the original purchaser of this product that it is free from defects in materials and workmanship for a period of 1 year from date of purchase.

If the unit becomes defective within the 1 year period, return the unit with original purchase receipt as your proof-of-purchase date. For return address information and instructions, go to <http://www.violight.com>.

Within 1 year of purchase, Violight Inc will repair or replace any part* of the product we deem defective at no charge, provided the product has not been abused, misused, or damaged, and was used according to these instructions.

*Defective bulbs will be replaced within 30 days of purchase.

This warranty gives you specific legal rights and you may also have other rights, which vary from state to state.

©2008 VIOlight Inc 888-882-2367 - customerservice@violight.com www.violight.com