

Limited 90 Day Warranty

VIOLight Inc. warrants to the original purchaser of this product that it is free from defects in materials and workmanship for a period of 90 days from date of purchase.

If the unit becomes defective within the 90 day period, return the unit along with the original packaging and purchase receipt as your proof-of-purchase date.

Please visit www.violight.com or call 888-882-2367 for return address and customer service information.

Within 90 days of purchase, VIOLight will repair or replace any part of the product we deem defective at no charge, provided the product has not been abused, misused, or damaged and was used according to the instructions.

This warranty gives you specific legal rights and you may also have other rights which vary from state to state.

vio^{light} duo



toothbrush sanitizer

Sanitizing is so easy!

Congratulations on your purchase of the VIOLight DUO Personal Toothbrush Sanitizer.

As a health-conscious consumer, you're aware of the benefits of keeping a clean toothbrush, and the importance of sanitizing it every day.

With your new VIOLight DUO Personal Toothbrush Sanitizer, you will be able to reduce your exposure to harmful germs and bacteria at home or while traveling.

VIOLight's germicidal ultraviolet light will eliminate up to 99% of the germs it comes into contact with on your toothbrush. You brush, rinse, and insert your toothbrush into the VIOLight – close the lid and sanitizing begins automatically.

Best of all, your toothbrush is stored in an enclosed, sanitized environment ready for the next time you brush.

Special offer for Violight UV Items!
FREE SHIPPING at violight.com when you enter coupon code "SHIP" at checkout.



Violight Personal/Travel Toothbrush Sanitizer.

A toothbrush case with built in UV power. Free Toothbrush included! \$29.95



Violight Countertop Toothbrush Sanitizer.

Sanitize 4 manual or electric toothbrushes! \$49.95



Violight Dental Spa.

UV + Sonic cleaning for dentures, retainers, mouthguards & more! \$59.95



Violight UV Wand.

Disinfect any surface with a wave of the wand. Great for countertops, keyboards, phones, & toys! \$59.95

Using Your VIOLight

Your VIOLight DUO requires 2 "AAA" batteries. The battery compartment is located on the back of your sanitizer.

1. Inserting Toothbrushes




Open the cover of your VIOLight DUO. Insert toothbrushes into rubber holders with the bristles facing the bulb.

2. Activating UV




Close cover to activate UV light. Your DUO shuts off automatically after 6 minutes – your toothbrush is fully sanitized.

3. Optional Wall Mount

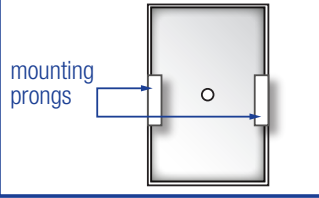

Affix the optional mounting bracket to your wall with the included adhesive pad, or with a screw (not included).

4. Mount & Remove

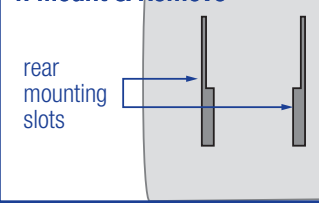

Align the mounting slots on the back of your VIOLight DUO with the prongs on the mounting bracket. Place Duo over the prongs, then slide down to lock onto bracket. Reverse steps to remove.

Safeguards & Precautions

1. Alkaline batteries are recommended for optimum performance.
2. Don't mix old and new batteries.
3. Do not immerse unit in water or liquid of any kind.
4. Avoid looking directly into the illuminated light when unit is turned on.
5. Parts of this device could present a choking hazard to young children. Keep out of reach of small children.
6. Unit is equipped with a safety mechanism that will not allow unit to go on when lid is open.
7. If the bulb should stay lit when lid is open, do not look directly into the lit bulb. Remove batteries, and visit www.violight.com for customer service assistance.

Cleaning & Care

1. Use a damp cloth to clean the inside and outside surfaces of your VIOLight Duo
2. The bulb in your VIOLight DUO is rated to remain effective for several years of use. If the bulb does not come on, simply replace the batteries in your DUO.

For questions, troubleshooting, or to order any of our UV products please call or email:

customerservice@violight.com

1-888-882-2367

www.violight.com

Copyright ©2008 Violight Inc.
VI0290-052808